

PCRmax Ltd TM qPCR test

shiga toxin (stx2b) producing Escherichia coli

E.coli generic detection of beta-
D-glucuronidase (uidA) gene
&
shiga toxin (stx2B) gene

150 tests

PCRmax 

For general laboratory and research use only

PCRmax 

Quantification of shiga toxin (stx2b) producing Escherichia coli genomes.

Advanced kit handbook HB10.07.08

Published Date: 20/09/2016

Introduction to shiga toxin (stx2b) producing *Escherichia coli*

Escherichia coli are one of many species of bacteria living in the lower intestines of mammals, known as gut flora. When located in the large intestine, it assists with waste processing, vitamin K production, and food absorption. Discovered in 1885 by Theodor Escherich, a German pediatrician and bacteriologist, *E. coli* are abundant: the number of individual *E. coli* bacteria in the faeces that a human defecates in one day averages between 100 billion and 10 trillion. However, the bacteria are not confined to the environment, and specimens have also been located, for example, on the edge of hot springs. The bacteria are Gram-negative, rod-shaped, flagellated and non-spore forming. Most strains are non-pathogenic but some cause food poisoning in humans with transmission largely being through the faecal-oral route. *E. coli* have a circular, DNA genome of approximately 4.6 Mb but also carry plasmids.

Shiga toxin-producing *E. coli* (STEC) are a form of enterohaemorrhagic *E. coli* that cause illness ranging from mild intestinal disease to severe kidney disease. The shiga toxin can cause haemorrhagic colitis, the source of the bloody diarrhoea associated with *E. coli* O157:H7 infections, as well as being responsible for haemolytic uremic syndrome (HUS). Shiga toxins derive their name from the organism where they were first classified, *Shigella dysenteriae*. When the shiga toxin is released, it can translocate to organs other than the digestive tract such as the kidneys and central nervous system. The ability of the shiga toxins to pass through cell barriers is possibly due to the increased permeability of the intestinal epithelial cells resulting from effects of the body's own immune system. The body increases permeability of cell barriers so that important cells of the immune system (neutrophils/PMN's) can reach the *E. coli* infection. Shiga toxin may use this opportunity to break through the walls of the digestive tract, enter the blood stream, and bind white blood cells for transport to locations such as the kidney or brain.

Enterohaemorrhagic *E. coli* are found in humans, cattle, and goats. There are a number of *E. coli* serogroups that produce shiga toxin such as O157:H7, O26, O111, and O103. Typical symptoms include severe abdominal cramping, sudden onset of watery diarrhoea, frequently bloody, and sometimes vomiting and a low-grade fever. Most often the illness is mild and self-limited generally lasting 1-3 days. However, serious complications such as haemorrhagic colitis, haemolytic uremic syndrome (HUS), or post diarrhoeal thrombotic thrombocytopenic purpura (TTP) can occur in up to 10% of cases. The incubation period ranges from 1 to 8 days, and transmission is predominantly through consumption of contaminated foods.

Specificity

The PCRmax qPCR Kit for shiga toxin (stx2b) producing Escherichia coli (E.coli_stx2b) genomes is designed for the in vitro quantification of E.coli_stx2b genomes. The kit is designed to have the broadest detection profile possible whilst remaining specific to the E.coli_stx2b genome.

The primers and probe sequences in this kit have 100% homology with a broad range of E.coli_stx2b sequences based on a comprehensive bioinformatics analysis.

If you require further information, or have a specific question about the detection profile of this kit then please send an e.mail to help@pcrmax.com and our bioinformatics team will answer your question.

Kit Contents

- **beta-D-glucuronidase primer/probe mix (150 reactions BROWN)**
FAM labelled
- **Shiga toxin 2B subunit primer/probe mix (150 reactions BROWN)**
FAM labelled
- **beta-D-glucuronidase positive control template (for Standard curve RED)**
- **Shiga toxin 2B subunit positive control template (for Standard curve RED)**
- **Internal extraction control primer/probe mix (150 reactions BROWN)**
VIC labelled as standard
- **Internal extraction control DNA (150 reactions BLUE)**
- **Endogenous control primer/probe mix (150 reactions BROWN)**
FAM labelled
- **RNAse/DNAse free water (WHITE)**
for resuspension of primer/probe mixes and internal extraction control DNA
- **Template preparation buffer (YELLOW)**
for resuspension of positive control templates and standard curve preparation

Reagents and equipment to be supplied by the user

Real-Time PCR Instrument

DNA extraction kit

This kit is designed to work well with all processes that yield high quality DNA with minimal PCR inhibitors.

Lyophilised 2x qPCR Mastermix

This kit is designed to be compatible with all commercially available Mastermixes that run with standard cycling conditions.

Pipettors and Tips

Vortex and centrifuge

Thin walled 1.5 ml PCR reaction tubes

Kit storage and stability

This kit is stable at room temperature but should be stored at -20°C on arrival. Once the lyophilised components have been resuspended they should not be exposed to temperatures above -20°C for longer than 30 minutes and unnecessary repeated freeze/thawing should be avoided. The kit is stable for six months from the date of resuspension under these circumstances.

If a standard curve dilution series is prepared this can be stored frozen for an extended period. If you see any degradation in this serial dilution a fresh standard curve can be prepared from the positive control.

PCRmax does not recommend using the kit after the expiry date stated on the pack.

Suitable sample material

All kinds of sample material suited for PCR amplification can be used. Please ensure the samples are suitable in terms of purity, concentration, and DNA integrity (An internal PCR control is supplied to test for non specific PCR inhibitors). Always run at least one negative control with the samples. To prepare a negative-control, replace the template DNA sample with RNase/DNase free water.

Dynamic range of test

Under optimal PCR conditions PCRmax E.coli_stx2b detection kits have very high priming efficiencies of >95% and can detect less than 100 copies of target template.

Notices and disclaimers

This product is developed, designed and sold for research purposes only. It is not intended for human diagnostic or drug purposes or to be administered to humans unless clearly expressed for that purpose by the Food and Drug Administration in the USA or the appropriate regulatory authorities in the country of use. During the warranty period PCRmax detection kits allow precise and reproducible data recovery combined with excellent sensitivity. For data obtained by violation to the general GLP guidelines and the manufacturer's recommendations the right to claim under guarantee is expired. PCR is a proprietary technology covered by several US and foreign patents. These patents are owned by Roche Molecular Systems Inc. and have been sub-licensed by PE Corporation in certain fields. Depending on your specific application you may need a license from Roche or PE to practice PCR. Additional information on purchasing licenses to practice the PCR process may be obtained by contacting the Director of Licensing at Roche Molecular Systems, 1145 Atlantic Avenue, Alameda, CA 94501 or Applied Biosystems business group of the Applied Biosystems Corporation, 850 Lincoln Centre Drive, Foster City, CA 94404. In addition, the 5' nuclease assay and other homogeneous amplification methods used in connection with the PCR process may be covered by U.S. Patents 5,210,015 and 5,487,972, owned by Roche Molecular Systems, Inc. and by U.S. Patent 5,538,848, owned by The Perkin-Elmer Corporation.

Trademarks

PCRmax™ is a trademark of PCRmax Ltd.

The PCR process is covered by US Patents 4,683,195, and 4,683,202 and foreign equivalents owned by Hoffmann-La Roche AG. BI, ABI PRISM®, GeneAmp® and MicroAmp® are registered trademarks of the Applied Biosystems (Applied Biosystems Corporation). BIOMEK® is a registered trademark of Beckman Instruments, Inc.; iCycler™ is a registered trademark of Bio-Rad Laboratories, Rotor-Gene is a trademark of Corbett Research. LightCycler™ is a registered trademark of the Idaho Technology Inc. GeneAmp®, TaqMan® and AmpliTaqGold® are registered trademarks of Roche Molecular Systems, Inc., The purchase of the PCRmax reagents cannot be construed as an authorization or implicit license to practice PCR under any patents held by Hoffmann-La Roche Inc.

Principles of the test

The kit contains two primer and probe sets. The uidA primer and probe set are designed to detect all E.coli sequences regardless of any other pathogenic markers that may be carried by the strain.

The stx2B primer and probe set is specific to the mobile genetic element that contains the Shiga toxin operon. Samples that test positive for uidA and stx2b and confirmed to be Shiga toxin-producing Escherichia coli (STEC). Samples that test positive for uidA but are negative for stx2B indicate that the sample contained an E.coli strain but not one that encodes the shiga toxin operon.

Real-time PCR

uidA and stx2B specific primer and probe mix is provided and this can be detected through the FAM channel.

The primer and probe mixes provided exploit the so-called TaqMan® principle. During PCR amplification, forward and reverse primers hybridize to the E.coli_stx2b DNA. A fluorogenic probe is included in the same reaction mixture which consists of a DNA probe labeled with a 5'-dye and a 3'-quencher. During PCR amplification, the probe is cleaved and the reporter dye and quencher are separated. The resulting increase in fluorescence can be detected on a range of real-time PCR platforms.

Positive control

For copy number determination and as a positive control for the PCR set up, the kit contains a positive control template. This can be used to generate a standard curve of uidA or stx2B copy number / Cq value. Alternatively the positive control can be used at a single dilution where full quantitative analysis of the samples is not required. Each time the kit is used, at least one positive control reaction must be included in the run. A positive result indicates that the primers and probes for detecting the target E.coli_stx2b gene worked properly in that particular experimental scenario. If a negative result is obtained the test results are invalid and must be repeated. Care should be taken to ensure that the positive control does not contaminate any other kit component which would lead to false-positive results. This can be achieved by handling this component in a Post PCR environment. Care should also be taken to avoid cross-contamination of other samples when adding the positive control to the run. This can be avoided by sealing all other samples and negative controls before pipetting the positive control into the positive control well.

Negative control

To validate any positive findings a negative control reaction should be included every time the kit is used. For this reaction the RNase/DNase free water should be used instead of template. A negative result indicates that the reagents have not become contaminated while setting up the run.

E.coli_stx2b DNA is known to be highly prevalent within the air and environment generally and the negative control may therefore give a late positive signal due to environmental contamination. The interpretation of results section of this handbook gives guidance on how to interpret results where environmental contamination is evident.

Internal DNA extraction control

When performing DNA extraction, it is often advantageous to have an exogenous source of DNA template that is spiked into the lysis buffer. This control DNA is then co-purified with the sample DNA and can be detected as a positive control for the extraction process. Successful co-purification and real-time PCR for the control DNA also indicates that PCR inhibitors are not present at a high concentration.

A separate primer and probe mix are supplied with this kit to detect the exogenous DNA using real-time PCR. The primers are present at PCR limiting concentrations which allows multiplexing with the target sequence primers. Amplification of the control DNA does not interfere with detection of the *uidA* or *stx2B* target DNA even when present at low copy number. The Internal control is detected through the VIC channel and gives a Cq value of 28+/-3.

Endogenous control

To confirm extraction of a valid biological template, a primer and probe mix is included to detect an endogenous gene. Detection of the endogenous control is through the FAM channel and it is NOT therefore possible to perform a multiplex with the *uidA* or *stx2B* primers. A poor endogenous control signal may indicate that the sample did not contain sufficient biological material.

Carry-over prevention using UNG (optional)

Carry over contamination between PCR reactions can be prevented by including uracil-N-glycosylase (UNG) in the reaction mix. Some commercial mastermix preparations contain UNG or alternatively it can be added as a separate component. UNG can only prevent carry over from PCR reactions that include deoxyuridine triphosphate (dUTP) in the original PCR reaction. PCRmax recommend the application of 0.2U UNG per assay with a 15 minute incubation step at 37°C prior to amplification. The heat-labile UNG is then inactivated during the Taq polymerase activation step.

Reconstitution Protocol

To minimize the risk of contamination with foreign DNA, we recommend that all pipetting be performed in a PCR clean environment. Ideally this would be a designated PCR lab or PCR cabinet. Filter tips are recommended for all pipetting steps.

- 1. Pulse-spin each tube in a centrifuge before opening.**
This will ensure lyophilised primer and probe mix is in the base of the tube and is not spilt upon opening the tube.
- 2. Reconstitute the kit components in the RNase/DNase free water supplied, according to the table below:**
To ensure complete resuspension, vortex each tube thoroughly.

Component - resuspend in water	Volume
Pre-PCR pack	
uidA primer/probe mix (BROWN)	165 µl
stx2B primer/probe mix (BROWN)	165 µl
Internal extraction control primer/probe mix (BROWN)	165 µl
Endogenous control primer/probe mix (BROWN)	165 µl
Pre-PCR heat-sealed foil	
Internal extraction control DNA (BLUE)	600 µl

- 3. Reconstitute the positive control templates in the template preparation buffer supplied, according to the table below:**
To ensure complete resuspension, vortex each tube thoroughly.

Component - resuspend in template preparation buffer	Volume
Post-PCR heat-sealed foil	
uidA Positive Control Template (RED) *	500 µl
stx2B Positive Control Template (RED) *	500 µl

* This component contains high copy number template and is a VERY significant contamination risk. It must be opened and handled in a separate laboratory environment, away from the other components.

DNA extraction

The internal extraction control DNA can be added either to the DNA lysis/extraction buffer or to the DNA sample once it has been resuspended in lysis buffer.

DO NOT add the internal extraction control DNA directly to the unprocessed biological sample as this will lead to degradation and a loss in signal.

- 1. Add 4µl of the Internal extraction control DNA (BLUE) to each sample in DNA lysis/extraction buffer per sample.**
- 2. Complete DNA extraction according to the manufacturers protocols.**

Real-time PCR detection protocol

1. **For each DNA sample prepare a reaction mix according to the table below:**
Include sufficient reactions for positive and negative controls.

Component	Volume
Lyophilised 2x qPCR Mastermix	10 μ l
uidA or stx2B primer/probe mix (BROWN)	1 μ l
Internal extraction control primer/probe mix (BROWN)	1 μ l
RNAse/DNAse free water (WHITE)	3 μ l
Final Volume	15 μl

2. **For each DNA sample prepare an endogenous control reaction according to the table below (Optional):**
This control reaction will provide crucial information regarding the quality of the biological sample.

Component	Volume
Lyophilised 2x qPCR Mastermix	10 μ l
Endogenous control primer/probe mix (BROWN)	1 μ l
RNAse/DNAse free water (WHITE)	4 μ l
Final Volume	15 μl

3. **Pipette 15 μ l of each mix into individual wells according to your real-time PCR experimental plate set up.**
4. **Prepare sample DNA templates for each of your samples.**
5. **Pipette 5 μ l of DNA template into each well, according to your experimental plate set up.**
For negative control wells use 5 μ l of RNAse/DNAse free water. The final volume in each well is 20 μ l.

6. **If a standard curve is included for quantitative analysis prepare a reaction mix according to the table below:**

Component	Volume
Lyophilised 2x qPCR Mastermix	10 μ l
uidA or stx2B primer/probe mix (BROWN)	1 μ l
RNAse/DNAse free water (WHITE)	4 μ l
Final Volume	15 μl

7. Preparation of standard curve dilution series.

- 1) Pipette 90µl of template preparation buffer into 5 tubes and label 2-6
- 2) Pipette 10µl of Positive Control Template (RED) into tube 2
- 3) Vortex thoroughly
- 4) Change pipette tip and pipette 10µl from tube 2 into tube 3
- 5) Vortex thoroughly

Repeat steps 4 and 5 to complete the dilution series

Standard Curve	Copy Number
Tube 1 Positive control (RED)	2 x 10 ⁵ per µl
Tube 2	2 x 10 ⁴ per µl
Tube 3	2 x 10 ³ per µl
Tube 4	2 x 10 ² per µl
Tube 5	20 per µl
Tube 6	2 per µl

8. Pipette 5µl of standard template into each well for the standard curve according to your experimental plate set up.
The final volume in each well is 20µl.

Amplification Protocol

Amplification conditions using Lyophilised 2x qPCR Mastermix.

	Step	Time	Temp
	UNG treatment (if required) **	15 mins	37 °C
	Enzyme activation	2 mins	95 °C
Cycling x50	Denaturation	10 secs	95 °C
	DATA COLLECTION *	60 secs	60 °C

* Fluorogenic data should be collected during this step through the FAM and VIC channels



** Required if your Mastermix includes UNG to prevent PCR carryover contamination

Interpretation of Results

Target (FAM)	internal control (VIC)	Positive control	Negative control	Interpretation
≤ 30	+ / -	+	-	POSITIVE QUANTITATIVE RESULT calculate copy number
> 30	+	+	-	POSITIVE QUANTITATIVE RESULT calculate copy number
> 30	-	+	-	POSITIVE QUALITATIVE RESULT do not report copy number as this may be due to poor sample extraction
-	+	+	-	NEGATIVE RESULT
+ / -	+ / -	+	≤ 35	EXPERIMENT FAILED due to test contamination
+ / -	+ / -	+	> 35	*
-	-	+	-	SAMPLE PREPARATION FAILED
+ / -	+ / -	-	+ / -	EXPERIMENT FAILED

Positive control template (RED) is expected to amplify between Cq 16 and 23. Failure to satisfy this quality control criterion is a strong indication that the experiment has been compromised.

*Where the test sample is positive and the negative control is positive with a Cq > 35, the sample must be reinterpreted based on the relative signal strength of the two results:

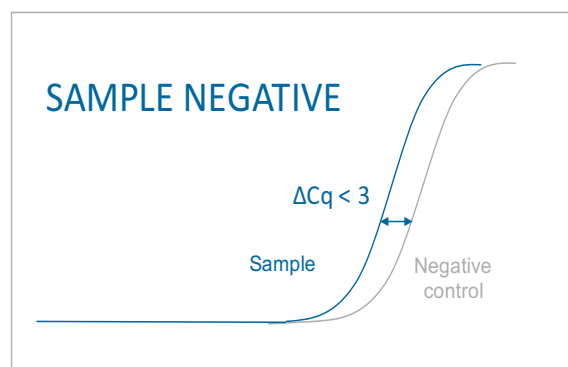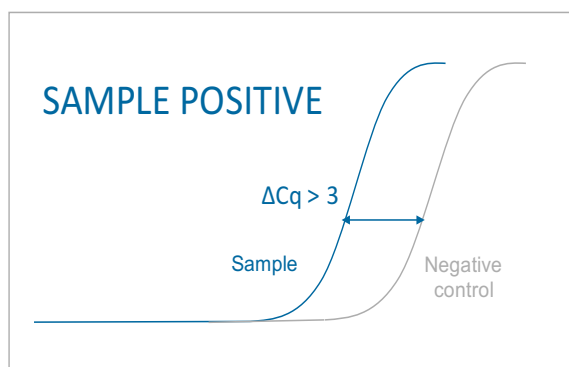

If the sample amplifies > 3 Cq earlier than the negative control then the sample should be reinterpreted (via the table above) with the negative control verified as negative.

If the sample amplifies < 3 Cq earlier than the negative control then the positive sample result is invalidated and a negative call is the correct result.

Internal PCR control

The Cq value obtained with the internal control will vary significantly depending on the extraction efficiency, the quantity of DNA added to the PCR reaction and the individual machine settings. Cq values of 28 ± 3 are within the normal range. When amplifying a *E. coli*_stx2b sample with a high genome copy number, the internal extraction control may not produce an amplification plot. This does not invalidate the test and should be interpreted as a positive experimental result.

Endogenous control

The signal obtained from the endogenous control primer and probe set will vary according to the amount of biological material present in a given sample. An early signal indicates the presence of a good yield of biological material. A late signal suggests that little biological material is present in the sample.